


Non-absorbable Surgical Contact Gripping Mesh
CONTACTO™ MESH
Sterilised Surgical Polyester / Polyurethane Mesh

DESCRIPTION


- Tridimensional Light “honeycomb” knitted structure with “grips”
- Polyester multi filaments impregnated with polyurethane (PUR)

Weight: 126 g/m² - Thickness : 1.09 mm
 Mean bursting strength: 314 kPa*
 Pore size: 1.48 x 1.32 mm - Porosity: 61 %
 Mean Tensile strength (Warp/Weft): 147/166 N
 Mean Elongation (Warp/Weft): 41/76 %
 *Max. intra abdominal pressure 24 kPa

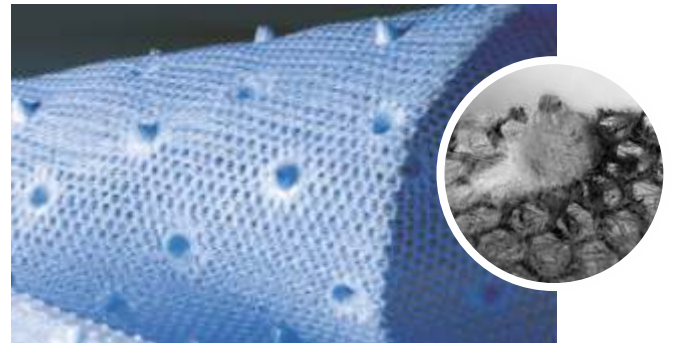

Advantages

- Self gripping prosthesis thanks to atraumatic “grips” on both sides
- Easy repositioning during surgery
- Porous structure to favour fast tissue ingrowth and colonization
- Shape memory for perfect use for laparoscopic or open surgery
- Elasticity and excellent multidirectional mechanical properties

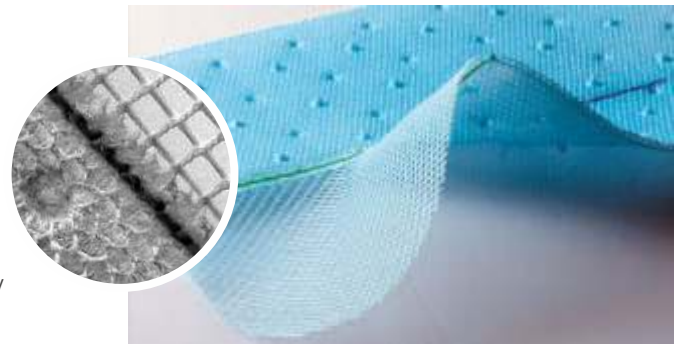

Indications

Laparoscopic Inguinal Hernias

Code	Size in cm	Shapes
CGM 1115 R	11 X 15	
CGM 1115 L	11 X 15	
CGM 1216 R	12 x 16	
CGM 1216 L	12 x 16	
CGM 1515	15 x 15	